

buy smart

ORDER FAX +49 2173 39800 20

Quantity	Art. No.	Product	Quantity	Art. No.	Product
MATRIX SYSTEM FOR ALLIGATOR ABUTMENTS					
<input type="checkbox"/>	2010.101	Equipment box with 3 tools (excluding consumables)	<input type="checkbox"/>	2010.731	Demounting tool for mounting inserts + model analogue reposition aid, blue
<input type="checkbox"/>	2010.601	Processing package Titanium 2 pcs. Titanium matrix housing incl. PEEK mounting insert black Retention inserts white, yellow, green Mounting collar silicone	<input type="checkbox"/>	2010.741	Mounting and demounting tool for retention inserts, brown
<input type="checkbox"/>	2010.611	Processing package PEEK 2 pcs. PEEK matrix housing incl. PEEK mounting insert black Retention inserts white, yellow, green Mounting collar silicone	<input type="checkbox"/>	2010.751	Matrix housing extractor gray
<input type="checkbox"/>	2010.710	Retention insert red PEEK, pull-off weight 4 pcs. extra-light	<input type="checkbox"/>	2010.703	Titanium matrix housing with attachment option incl. PEEK mounting insert 4 pcs.
<input type="checkbox"/>	2010.711	Retention insert white PEEK, pull-off weight 4 pcs. light	<input type="checkbox"/>	2010.701	Titanium matrix housing incl. POM mounting insert 4 pcs.
<input type="checkbox"/>	2010.712	Retention insert yellow PEEK, pull-off weight 4 pcs. medium	<input type="checkbox"/>	2010.702	PEEK matrix housing incl. POM mounting insert 4 pcs.
<input type="checkbox"/>	2010.713	Retention insert green PEEK, pull-off weight 4 pcs. strong	<input type="checkbox"/>	2010.725	POM mounting insert white 4 pcs.
<input type="checkbox"/>	2010.714	Retention insert blue PEEK, pull-off weight 4 pcs. extra-strong	<input type="checkbox"/>	2010.721	Model analog Aluminium 4 pcs.
<input type="checkbox"/>	2010.715	Retention insert black PEEK, pull-off weight 4 pcs. ultra-strong	<input type="checkbox"/>	2010.722	Forming / fixing matrix, red PEEK 4 pcs.
			<input type="checkbox"/>	2010.723	Processing spacer, white POM C 4 pcs.
			<input type="checkbox"/>	2010.724	Mounting collar silicone 10 pcs.

Art. No. 90533 - DRS International GmbH - ISSUE December 2016



Novaloc™
s w i s s m a d e

MATRIX SYSTEM

The **Novaloc™** matrix system with its newly developed technology is a prefabricated connector for fixing removable dental prostheses on ALLIGATOR abutments.



Original Novaloc™ available by DENTAL RATIO NEW



DRS International GmbH
Raiffeisenstr. 15 a
40764 Langenfeld
Germany
Phone +49 2173 39800 0 Fax +49 2173 39800 20

Standard Express

Customer ID
Clinic/Lab
Name
Street
Zip Code/City
Phone
E-Mail

www.dental-ratio.de · www.valoc.ch



DRS International GmbH
Raiffeisenstr. 15 a
40764 Langenfeld
Germany
info@dental-ratio.de
www.dental-ratio.de
Phone +49 2173 39800 0
Fax +49 2173 39800 20


10 reasons why you choose Novaloc™

1 Matrix housing made of titanium or PEEK



Matrix housing in beige PEEK; the solution when there is very little space or when not using metal is a prerequisite

2 Stress-free mounting and removing of retention inserts

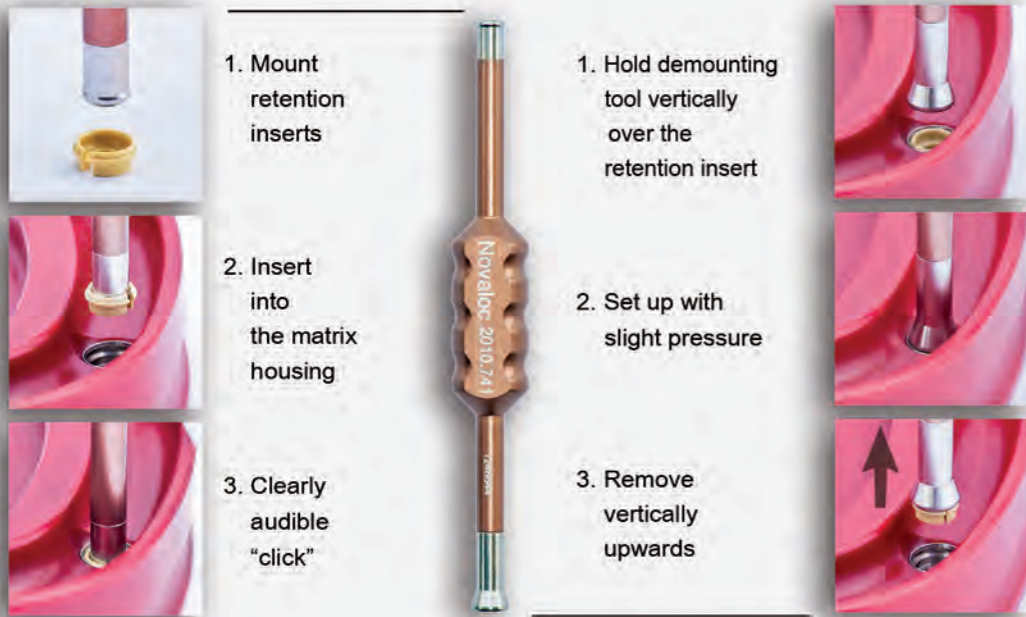

1. Mount retention inserts

2. Insert into the matrix housing


3. Clearly audible "click"

1. Hold demounting tool vertically over the retention insert

2. Set up with slight pressure

3. Remove vertically upwards

3 The ingenious Novaloc™ model analogue reposition aid

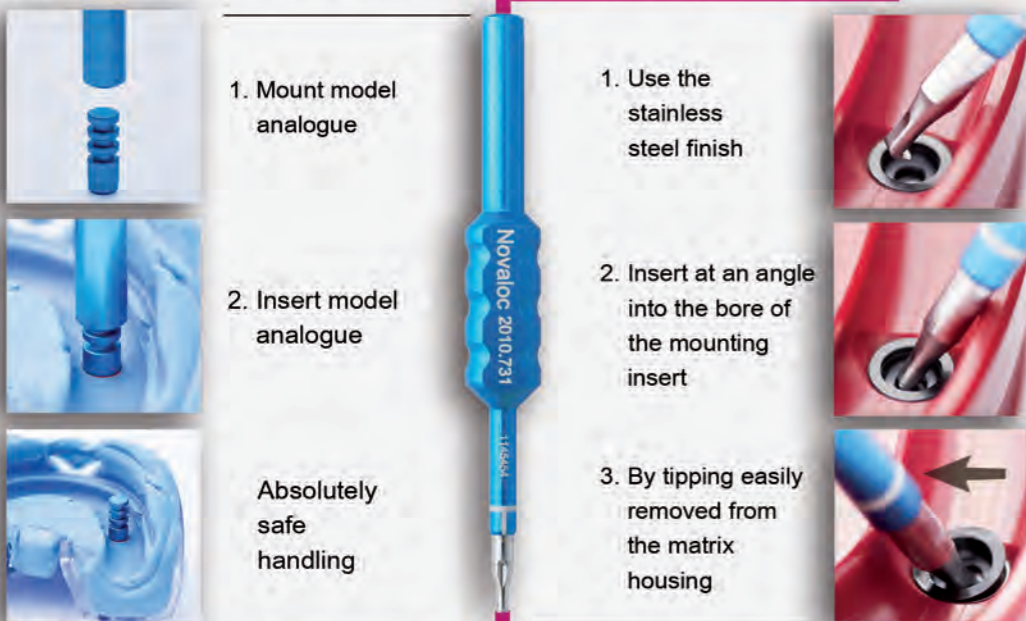

1. Mount model analogue

2. Insert model analogue

Absolutely safe handling

1. Use the stainless steel finish

2. Insert at an angle into the bore of the mounting insert

3. By tipping easily removed from the matrix housing

4 The clever demounting tool for mounting inserts

5 The user-friendly Novaloc™ equipment box



- No caps in the inventory
- Ensures overview and order
- No mixing when the cover is closed
- Ordering number legible in the cover
- Three special instruments for all sorts of manipulations

6 The "5-seconds" matrix housing extractor



1. Heat extraction head

2. Heat matrix housing briefly

3. Apply leverage to remove matrix housing

Special design: No more hot finger!

7 Novaloc™ only requires 4 replacement males



Retention value **light**

Retention value **medium**

Retention value **strong**

Retention value **extra-strong**

Easy to understand colour code according to retention value

8 New forming/fixing matrix for a safe and accurate positioning in the mouth



Now clearly smaller and more precise than those forming matrices commercially available until now

Multivalent can also be used as a fixing matrix for templates, bite registers and rails



9 Processing spacer now with blocking characteristics



A 5° inclination for the shaping of a self-retaining box for matrix attachment



For additional mechanical blocking

10 Clearly less wear and tear as well as servicing expense at the retention inserts

- Novaloc™ retention inserts need no central retention element (considerable advantage with regard to damage)
- PEEK has better physical properties and is considerably more hygienic than Nylon



Conventional damage to Nylon inserts commercially available until now